

Robotic Endoscope Holder MTG-H100





Robotic Endoscope Holder

MTG-H100



Wymiary: 1162(L) x 270(W) x 810(H) mm
Waga: Robot: 19.5 kg
Wózek: 92.5 kg
Zasilanie: AC110-230 V, 50 Hz/60 Hz, 2.5A
Akcesoria: jednorazowe zapięcia (MAS0001)
wózek(MAS0002)
sterylne pokrowce na ramię (MAS0003)
* Żaden element urządzenia nie może być
sterylizowany

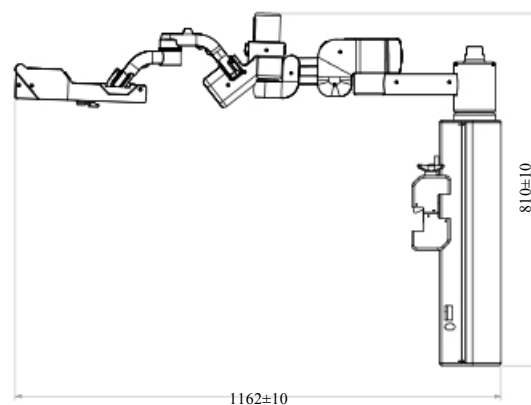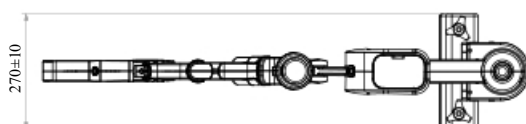



Kontrola nad kamerą i optyką

Robotyczne ramię (MTG-H100) służy do zamocowania kamery wraz z optyką by stabilizować obraz w trakcie zabiegu laparoskopowego. Konstrukcja urządzenia pozwala na montowanie ramienia bezpośrednio do stołu operacyjnego lub na specjalnym wózku.



Stabilizacja obrazu

MTG-H100 oferuje precyzyjne pozycjonowanie kamery, by zapewnić płynność obrazu. Eliminuje koszty związane z asystą do zabiegu oraz problemy z komunikacją.

Unikalny moduł RCM

The Remote Center of Motion module




Sześć zakresów kontroli

Aktywator nożny posiada sześć zakresów kontroli ruchu endoskopu w zakresach: przybliżanie/oddalanie, w górę/w dół, w prawo/w lewo. Porusza się z wyjątkową płynnością i precyzją.



Sterylna osłona

W trakcie zabiegu ramię jest chronione przez dedykowane osłony by zapobiec kontaminacji.

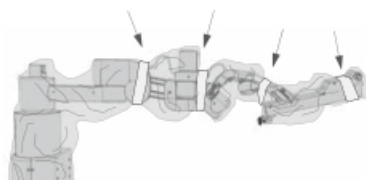

A Participant in Improving Quality of Human Life

Subsidiaries & R&D Centers

HIWIN GmbH
OFFENBURG, GERMANY
www.hiwin.de
www.hiwin.eu
info@hiwin.de

HIWIN Schweiz GmbH
JONA, SWITZERLAND
www.hiwin.ch
info@hiwin.ch

HIWIN KOREA
SUWON, KOREA
www.hiwin.kr
info@hiwin.kr

HIWIN JAPAN
KOBE · TOKYO · NAGOYA · NAGANO ·
TOHOKU · HOKURIKU · HIROSHIMA ·
KUMAMOTO · FUKUOKA, JAPAN
www.hiwin.co.jp
info@hiwin.co.jp

HIWIN FRANCE
ECHAUFFOUR, FRANCE
www.hiwin.fr
info@hiwin.fr

HIWIN CHINA
SUZHOU, CHINA
www.hiwin.cn
info@hiwin.cn

HIWIN USA
CHICAGO · SILICON VALLEY, U.S.A.
www.hiwin.com
info@hiwin.com

HIWIN s.r.o.
BRNO, CZECH REPUBLIC
www.hiwin.cz
info@hiwin.cz

Mega-Fabs Motion System, Ltd.
HAIFA, ISRAEL
www.mega-fabs.com
info@mega-fabs.com

HIWIN Srl
BRUGHERIO, ITALY
www.hiwin.it
info@hiwin.it

HIWIN SINGAPORE
SINGAPORE
www.hiwin.sg
info@hiwin.sg

HIWIN TECHNOLOGIES CORP.

No. 7, Jingke Road,
Taichung Precision Machinery Park,
Taichung 40852, Taiwan
Tel: +886-4-23594510
Fax: +886-4-23594420
www.hiwin.tw
business@hiwin.tw

- HIWIN is the registered trademark of HIWIN Technologies Corp.. For your protection; To avoid counterfeit products, be certain you are buying genuine HIWIN products before purchase.
- Actual products may be different from the specifications and photos in this catalog. The differences in appearances or specifications may be caused by, among other things, product improvements.
- HIWIN will not sell or export those techniques and products restricted under the "Foreign Trade Act" and relevant regulations. Any export of restricted products should be approved by competent authorities in accordance with relevant laws, and shall not be used to manufacture or develop nuclear, biochemical, missile and other military weapons.
- HIWIN website for patented product directory: http://www.hiwin.tw/Products/Products_patents.aspx